## NEATECH

Job Series















The JoB WALKER is a walker that allows people with temporary or permanent disabilities to access harsh environments such as the beach or snow. This is obtained thanks to the use of patented innovative solutions from Neatech, as the wheels made of plastic material with controlled deformation that do not sink on soft soils. The frame is made of light aluminum alloy, a material that guarantees a high resistance both to bending stress caused by the support of the user's weight and to the pushing strength, and to the saline environment corrosion. The JoB WALKER is a walker for people with physical disabilities or simply with mobility problems due to the age, who wish to deal with soft and irregular soils on whose it would be impossible to use traditional supports. The wheels, always ensuring the greatest possible surface for resting on the ground, allows ambulation not only of sand or snow, but also on other types of non-cohesive soils such as meadows, all these make this product a valid aid usable in all seasons and in all contexts. The folding seat allows the user to sit down for a stopover on the way, or to enjoy the sightseen. The back and seat are made of sturdy fabric characterized by rapid drying that does not incorporate the grains of the sand.

## **FEATURES**

resisting materials



inconsistent grounds


comfortable

## **TECHNICAL SPECIFICATIONS**

| Total length              | 380 mm               |
|---------------------------|----------------------|
| Total width               | 900 mm               |
| Minimum transport length  | 900 mm               |
| Seat height               | 500 mm               |
| Total weight              | 5,5 kg + 3kg x wheel |
| Maximum user weight       | 100 kg               |
| Seat width                | 460 mm               |
| Seat depth                | 390 mm               |
| Backrest width            | 460 mm               |
| Backrest height           | 445 mm               |
| Backrest angle-adjustable | fixed                |
| Legrest angle-adjustable  | fixed                |
| Warranty                  | 2 years              |
| CE EN12182                | yes                  |
|                           |                      |

maximum ease of use

## **AVAILABLE COLOURS**

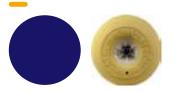

**NEATECH.it s.r.l.** 

fax e-mail info@neatech.it